



Nebulizer Set-up



- Twist off the top of one saline vial  and empty into chamber.
(one vial is equivalent to 5ml of 0.9% NaCl)
- Assemble the nebulizer by placing the baffle inside the chamber.
(Note: baffle may not be exactly as pictured to left)
- Screw on nebulizer cap.
- Attach mask to nebulizer cap.
- Attach tubing from bottom of chamber (air-inlet connector) to nebulizer machine outlet.
- Turn on compressor and place mask over mouth and nose

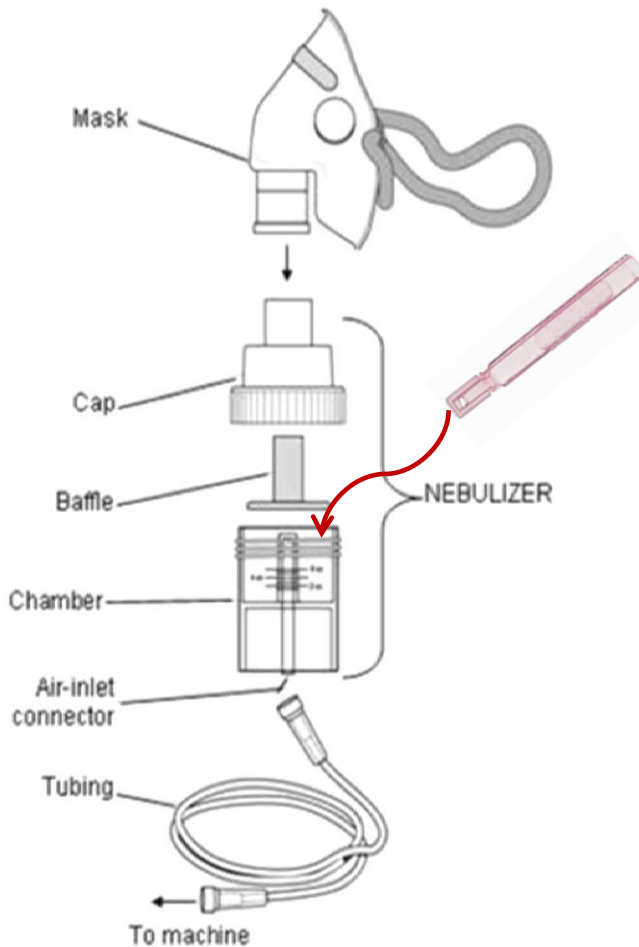

- Maintain upright position to prevent saline from spilling into mask
- Treatment takes approximately 10 minutes or when mist dissipates.

Compressors



- To clean: wipe mask after each use. Rinse the cap, baffle, and chamber to remove salt residue. Air dry.
(Wipe surface of compressor and tubing as needed.)

For more information:

www.technologyforliving.org